
Care lift






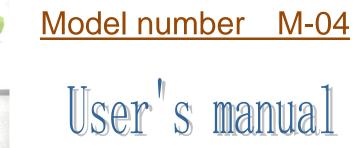


1
2
4
6
7
11
12
14
16


1 Accessories and Nomenclature

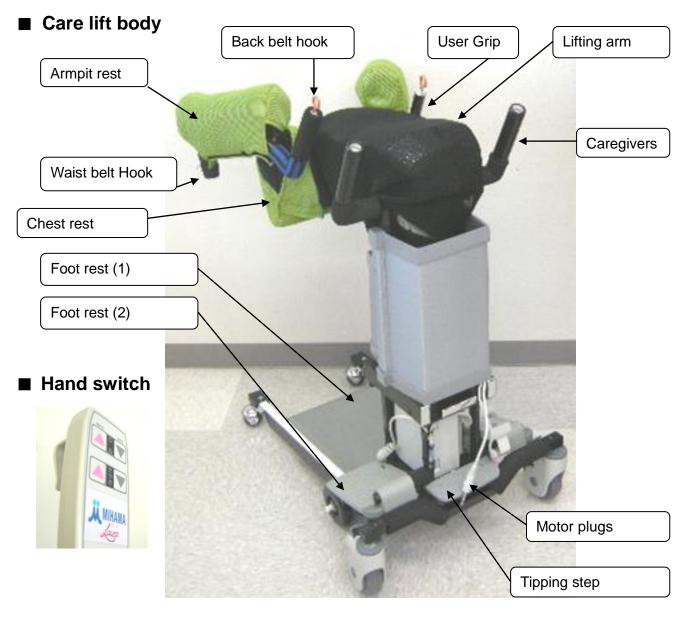

■ Two-piece unit

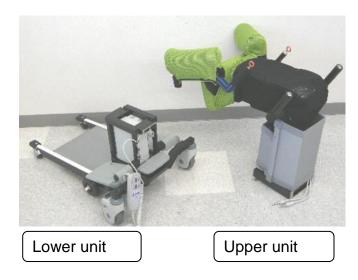


Accessories

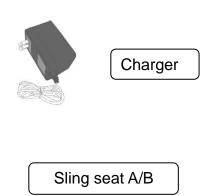

2 Assembling method


■ Assembly of the movement for the lift as "lower unit" of the "upper unit"

① Remove the "mounting bar" of the upper unit.

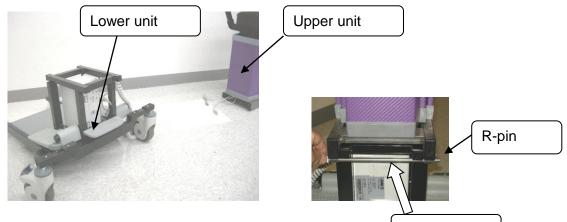



Mounting bar

② Lock the caster, and insert Tilt the upper frame diagonally.

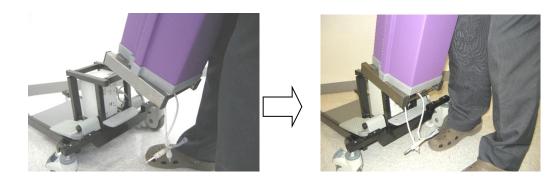

③ Insert the "mounting bar" of the upper frame, and secure it with the R pin. [Can be in either direction]

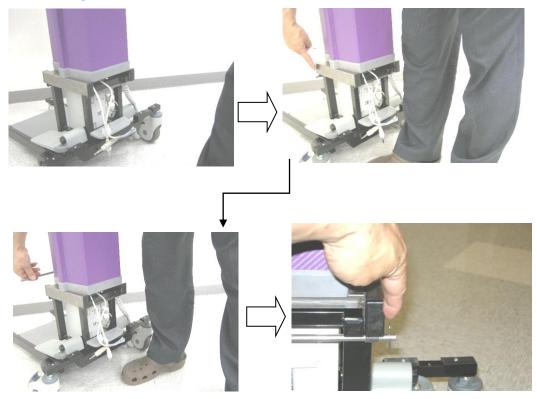

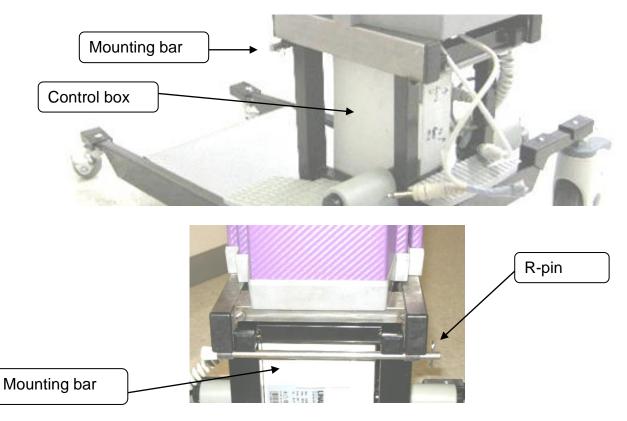


④ Make sure the upper frame and the lower frame has been fixed.


⑤ Connect the "Motor plug" of the upper frame to the "control box" of the lower frame.



"Insert the motor plug 1 (red mark) to the control BOX (red mark)."

"Insert the motor plug 2 (blue mark) to the control BOX (blue mark)."

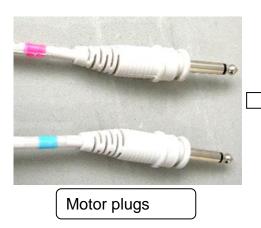

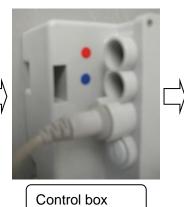

Connect the plug 1

3 Attention for Safety.

Because it is important information about safety, please be sure amulet.

(Description of the mark)


X the extent of the resulting harm or damage when it was the wrong use

are classified in the following display marks, it is described.

. Warning is "People are assuming is a possibility that death or serious injury" is the content.

Caution is "People only possibility or property damage injury is assumed likely to occur," it is the contents.

After you have finished reading, in the place where the person who is your seen at any time, please be sure to store.

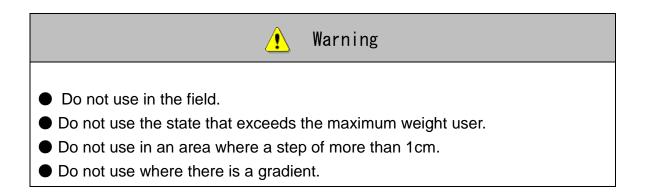

Do not stand up on top of the lifting arm and chest addressed.
• it is not running in a state where the lock casters. Moreover, it is not locked
during travel.
\bullet sudden acceleration, sudden stop, it is not a sudden change of direction.
● it is not allowed to use in small children.
ullet in a state in which the person has got, it does not lift the product.
• it is not or embarked the body during running.
ullet Do not push the body from the side while applying the brake lock lever.
Do not two-seater.
• speed to move the wheel to ensure a safe speed the step was also included.
sure you lock the casters when getting on and off.
ullet When you are traveling, or stepped foot on the wheel, and be careful not to get
caught or the hand.
 In steps and grooves, and enters vertically.
such as the dancer uneven, caregiver support the body.
Make sure not teetering hook.
• "breast addressed" and "side addressed to" Make sure that you are fixed with
screws.
Failure to comply with the above, injury, fall, it may cause a fall accident.
Keep away from the body to the fire.
Failure to comply with the above, it becomes body burning in a fire.
Do not touch the power plug with wet hands.
Do not charge in a place that is exposed to water.
■ If you do not adhere to the above, there is a risk of electric shock.
A Caution


- **L** Caution
- Do not apply heavy objects, such as handle armpit pad adjustment lever.
- disassembly, repair, is not carried out remodeling.
- Do not leave in direct sunlight.

■ Failure to comply with the above, the product of degradation is accelerated, it may cause malfunction or failure.

4 Each part operation method


1 Caution Make sure each part is working properly and fixed before

Hand switch

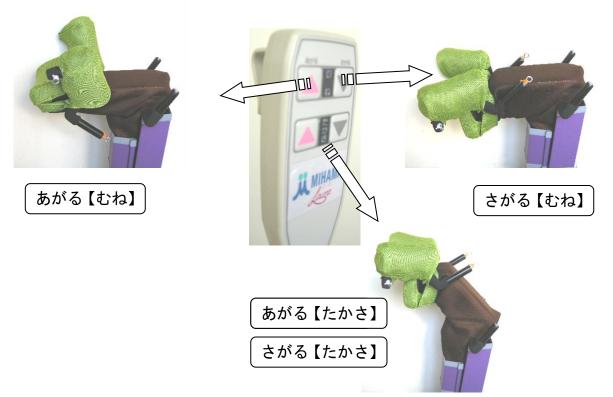

Overcoming step

• caregiver I floated the wheel stepped on tipping step.

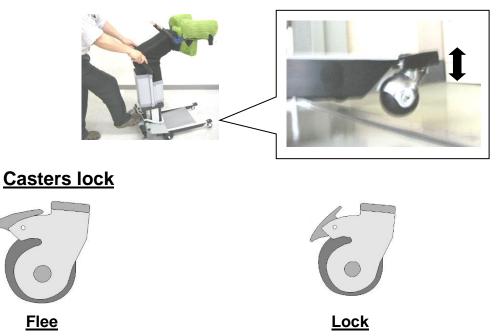


 \circ Lock the wheels, to secure the body.

 \circ left diagram in a free state, it is in the state that the right figure is locked.


5 Use procedure




-<

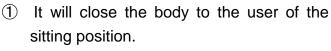

In this case, in accordance with the user, Armpit rest, Chest rest, the Lifting arm will keep adjusted with a Hand switch to fit the body.

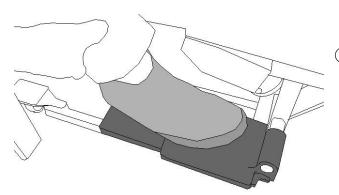

② Left and right feet, respectively it will put on Foot rests (1).

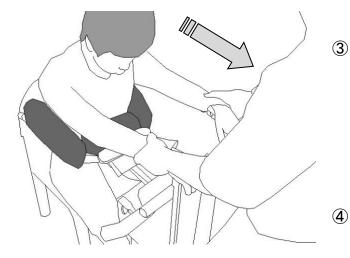

③ Place the both sides of the user Armpit rest, it is adhered to a position where the Chest is hit in the Chest rest. Then, lock the Casters to secure the body.




\rm L Caution

Once you lock the caster, please be sure to confirm or not to move the body.

Sling seat A


• Install the left and right of the snap (orange belt) on the back belt Hook by turning them back.

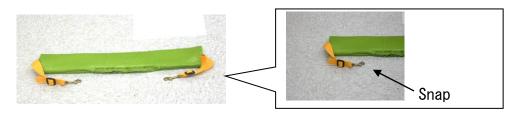

Sling seat B (Use from the chair)


• Install the left and right of the snap (orange belt) on the back belt Hook by turning them back.


• attached the left and right of the snap (white belt) to the waist belt Hook.

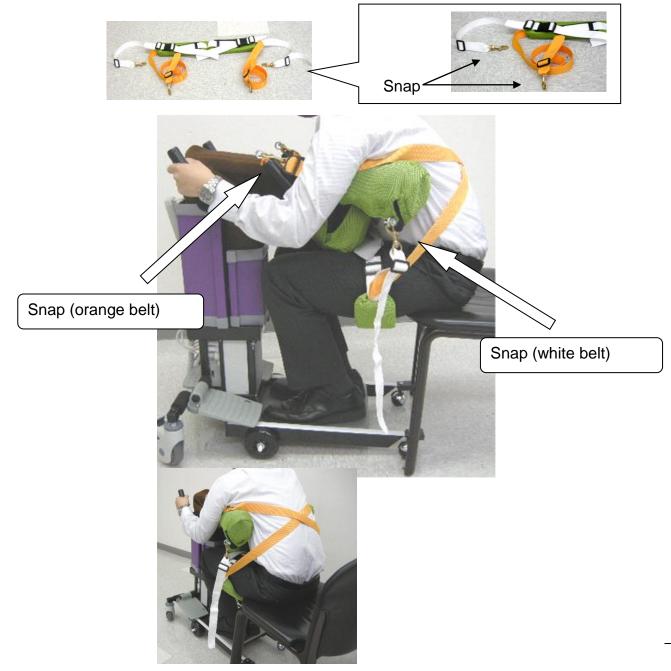

Sling seat B (Use from the floor)

- belt I attached so that there is no slack.
- it will put your feet on the footrest (2) at the position of the D.

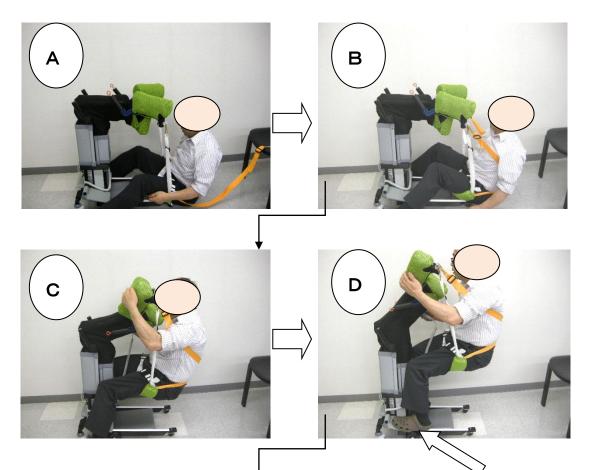

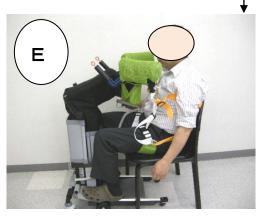

: Position of head

We will put your feet on the footrest (2)



If the height of the waist is low, caregiver it will raise the height by pulling slightly behind the hips.

🔥 Warning

Please be caregiver support the always body behind. Sling seat, please attach a position close to the ass.

1 Caution

• When operating the remote control switch, make sure there is no jamming. Also, if the user has complained of a malfunction, please discontinue use immediately.

If there is a danger that fall behind, use the "slipping prevention belt" attached.
 Sling sheet.

Please be mounted at a position close to the ass.



Please be mounted at a position close to the ass sling



Using the charger

 Please insert the plug of the charger to the charging port of the control BOX of lower unit

 \circ Please charge at 100V power supply for home.

 Sound of "Pihh-----" at the time of operation is small remaining amount of the charge comes from the control BOX.

• Please refrain from operation during

6 Precautions

■ Request to better to assistance In nursing lift is that I would like you to keep in mind when traveling with assistance. Please observe always. Please be sure to avoid the use in the following conditions.

- Steep slopes.

- Falling in fear and instability to become the location of the step.

- Such as the groove of a certain place getting on and off.

- The risk is considered location.

- Part of the person's body to sit while driving or touch (hands and feet) and the clothing on the body and the wheels, please check the operation so that it is not or pinched.

- If you want to assist towards the dementia, cushion, please pay close attention because it is conceivable to swallow such parts.

hope of slipping down to prevent belt wear From the care lift, it will prevent the full shift.
In such a reduction and joint contractors standing difficult people and the maintenance of muscle strength, people who do not try to stand up.

■ note of when you do not want to use the the care lift And long-term care it is that I would like you to note the location to put the lift. Please observe always. And outdoor and Direct sunlight (car included).

- Damp and were exposed to rain.

- The middle of a steep slope. Stove where you use the fire becomes hot, such as. And children a place that could be a prank.

- A high-traffic location.

Before you use

% Before using the care lift, please do the inspection of various parts for safety.

- Wheel confirmation of Whether there is a dent in the tires.

- The tire is no scratches, if there is any foreign matter. Confirmation of the lock mechanism with lock casters you can lock or release

- Or the wheel does not rotate in a state where a lock.

🔥 Warning

Please perform regular inspection and maintenance. The cause of the accident.

Periodically check the following things, please put in place. Upset in the case of a notice, please contact your dealer as soon as possible.

Department	Confirmation	Enabled	Check
position	item		
It does not move	The battery is	Please Charge.	
when you press the	low.		
remote control	There are beeps	Please Charge.	
switch.	that Sound of		
	"Pihh"		
	Life out of the	Please replace the battery.	
	battery.	% life standard of the battery is	
		about three years, it depends on the	
		method used.	
Caster of the lock is	There are foreign	Please remove.	
not working.	body in the tire.		
Body does not	Casters lock	Please unlock the caster.	
move.	state.		
The movement of	There is a	Foreign matter of pinching and each	
the body is	problem with	part of the flaws, please check the	
abnormal.	foreign objects or	loose bottles and nut portion.	
	component		
	mounting.		
Body is dirty.	There are dirt,	Then taking the dirt with a neutral	
	such as mud and	detergent, please wipe off the	
	dust.	detergent with a cloth and a dry	
		cloth containing clean water.	

Handling of battery

Battery and Notes on charging

The battery of the exchange, please contact us to like our company or dealer.

-When operating in charge in the remote control switch, battery life please be avoided will be shorter.


- In the code of charger wheel and heavy object, please do not step on the pinch.

-Remain connected to the code of the charger, please do not move the motorized ERROR device.

When you disconnect the power-plug, please do not pull out with a code.

8 Specifications

Product Name	Laura			
Type number	M-04			
The maximum user	90kg			
weight [Guideline	However, 70kg If you ride from the floor			
]	※ Please do not exceed this.			
Proper height [135 to 185 cm			
guide〗	Please do not exceed the height from 145 to less than 175			
	cm.			
Power supply	Rated voltage and	AC 100v · 50 / 60Hz		
	frequency	6VA within		
	Rated power	Charger (accessories)		
	consumption			
	Charging			
	apparatus			
Battery	Type formula	Two sealed lead-acid battery ×		
	Rating capacity	12vDC / 2.9Ah 2 pieces ×		
Motor	Type formula	Linear actuators: LA31		
	Rated voltage and	DC 24v · 2A		
	current	Duty ratio of 10%		
	Continuous use	To rest for 90 seconds when you use		
	time	※ 10 seconds		
	Lifting time	About 20 seconds		
Charging device	Input	AC 100v · 50 / 60Hz		
	Output	DC 26.5v · 200mA		
Minimum turning	45 cm			
radius [about]				
Minimum ground	3 cm			
height [about]				
Weight 〖about〗	28kg			
Material	Steel			
Paint	Black paint			
Caster	4 x diameter 10cm, 2 x diameter 5cm			
Production country	made in Japan			

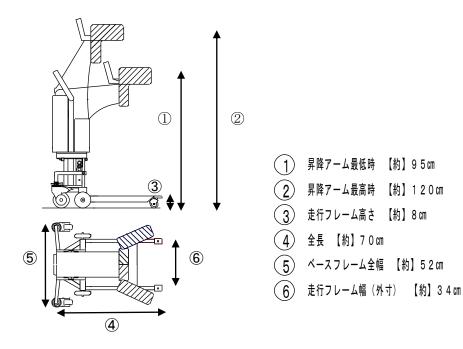

9 Warranty and after-sales service

When you receive the after-sales service or repair,

After-sales service, or when it is subject to repair, to "your address," "Your name", "phone number", "product name", "date of purchase", "purchased sale shop name," "failure or abnormal contents" to the contact below please let me know.

Mihama.co.ltd

Zip-code: 739-0046 Create core 1-4, Kagamiyama 3-13-60, Higashihiroshima-shi, Hiroshima-ken, JAPAN E-mail: info@mihama-care.jp

If you are past the warranty period

Repair, will be repaired by hope for a fee. However, the product may not be repaired.

- Overall maintenance (Not including battery)
- "Cost: about ¥ 40,000" ·· [lower unit] or [Upper unit]
- Please ask the amount of detail in the E-mail
- You may want to change the price and delivery time.



1 6